

Improve Stem Cell Yield

Sci-tive

Stem Cell Investigations - Total In Vitro Environment

Sci-tive has been designed to mimic in-vivo conditions providing a continuous cell culture environment which eliminates cellular stress linked to variations in temperature, pH and oxidation yielding better stem cell lines.

🌀 Complete “Lab in a Box”

Reduces need for laminar air flow cabinets, stereo microscopes, CO₂ incubators

🌀 Stress-Free Cell Manipulation

Now Stem Cell Biologists can work under stable O₂, CO₂, pH, temperature and humidity for derivation, manipulation and culturing steps

🌀 Direct-Hand Access into Workstation

Using sterilised and user friendly cuff & sleeve system (Ezee Sleeve™)

🌀 Imaging Ready


A digital stereo microscope or real-time live cell imaging system is easily accommodated

🌀 Aids Compliance with EU TCD (2004/23/EU)

VOC filtration and HEPA filtration create a biologically safe working environment

🌀 CE (1293) Approved (MDD 93/42/EEC)

Compliant with European EMC and electrical safety tests

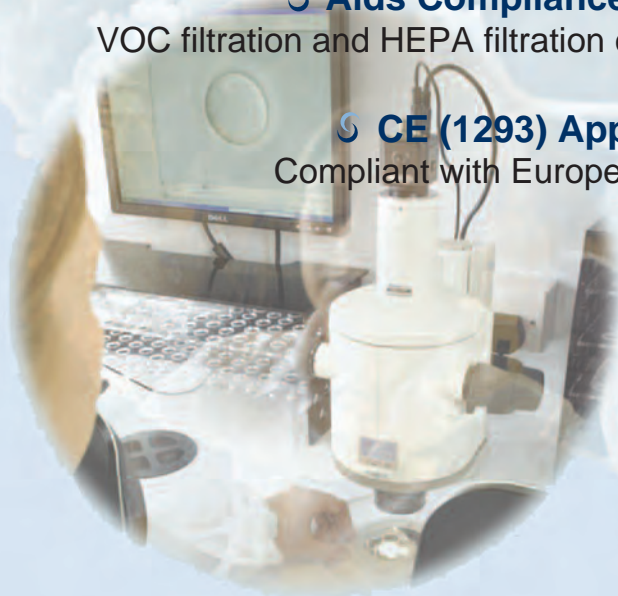

innovation in incubation

Product Specifications

External Workstation Dimensions

- Height of Unit (excluding stand) 1024mm
- Width of Unit (side to side) 1650mm
- Depth 760mm
(additional 150mm required behind workstation for piping)

Internal Workstation Dimensions

- Height 590mm
- Width 1174mm
- Depth 570mm

Workstation Weight

- Weight approx. 230kg

Interlock Dimensions

- Height 210mm
- Width 270mm
- Depth 200mm

Atmosphere Control

- Temperature Ambient + 5°C - 45°C
- Humidity Control Ambient - 85%
- O₂ Control 0.0% - 20.9%
- CO₂ Control 0.0% - 30.0%
- Integrated HEPA Filtration System
- Activated Carbon Filtration System

Alarms

- Temperature - Visual and Audible
- Humidity - Visual
- Gas Low Pressure - Visual and Audible

Standard Accessories

Workstation

- Stand
- 2 Racks with Perforated Shelves to Optimise Homogeneity
- Direct-hand Access using Ezee Sleeve™ and Ezee Pleez System
- Removable Front with UV Filtration Screen
- Waste Trays
- Gas Tight Waste Port

System PC and Software

- On-screen Keyboard
- Independent Touch Pad Mouse
- Windows Vista
- Flat Screen Monitor
- Imaging Camera
- Cronus3™ Imaging Software from RI

Gas Control Software

- USB Communications Port
- 7 Days Continuous Storage of Event Log Data
- One Touch Calibration (O₂)
- On-Screen Fault Assistance

Microscope

- Nikon AZ100 MZ Microscope fitted with triple nose piece and a single 2x magnification lens is supplied as standard which offers 46mm of working height. Other lenses can be purchased separately (choice is 1x or 4x magnification)

Sci-tive

Stem Cell Investigations - Total In Vitro Environment

The Ruskinn brand was founded in 1993 and rapidly became established as one of the world's leading suppliers and manufacturers of anaerobic and modified atmosphere workstations. Ruskinn's purpose built factory is located within the Sony Technology Centre near Pencoed in South Wales, UK. Ruskinn's range of high quality workstations is available worldwide through a global network of carefully selected distributors.



Digital, Stereo Microscope Options

- Nikon AZ100 MZ (as shown)
- Leica M 205 C
- Olympus

NEW Optional Accessories

- Remote Monitoring

For more information contact:

Ruskinn Technology Ltd
Suite 3 Technium Digital, Sony Technology Centre,
Pencoed CF35 5HZ, UK

Tel: +44 (0) 1656 868540 • Fax: +44 (0) 1656 868541

www.ruskinn.com • sales@ruskinn.com

Distributed by: